

CHAIR OP-S8

EXCELLENCE IN ENT CHAIRS

- Resistant
- Comfortable
- Totally electrical



OPTOMIC®

producing quality



producing quality

[OP-S8]

FOLDING ARMRESTS


They are composed of a stainless steel axle padded with washable polyurethane. They fold down completely to provide more comfort for patients with reduced mobility and for lateral observations and practices.



STRONG HEADREST

There are two models available, both of them perfectly designed to keep the head of the patient in the required position during the examination.

Adjustable in height, it incorporates a quick removal system which facilitates the execution of positional tests and routine practices.

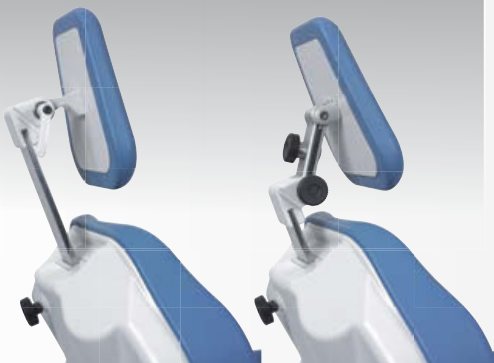

COMFORTABLE FOR DOCTOR AND PATIENT

AN ERGONOMIC PRODUCT

The OPTOMIC chairs design is based on the search for the comfort of the patient and ergonomics during exploration. This is the reason the geometry of its structure adapts to achieve the optimum distance between patient and doctor in every step of examination and treatment.



ROTATION AND BRAKE

ACCESSIBLE SYSTEM


The patient can be turned with great safety, stability and smoothness due to its high resistance bearings of conical rollers.

It has a double lever braking system, placed on either side of the seat, to avoid the doctor wasting time moving around the patient.

The 300° turn facilitates the examination of both ears.



CHAIR OP-S8

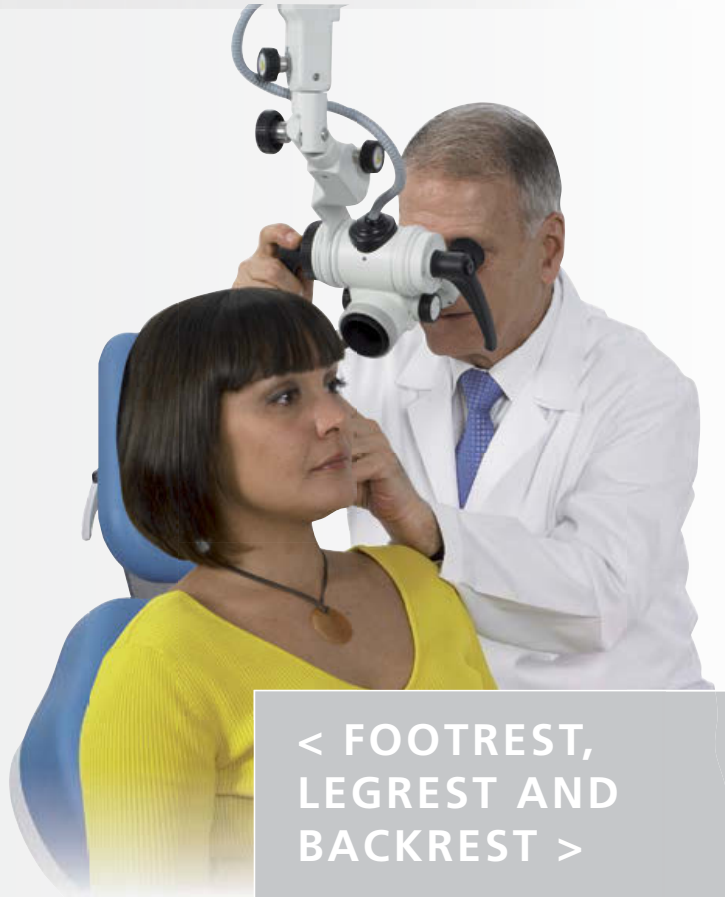

< FOOTREST,
LEGREST AND
BACKREST >

MOTORIZED ELEVATION SYSTEM

A PRECISE SYSTEM

A pedal facilitates the perfect control of the chair. It activates the up /down and inclination movements. It has two memories for user personalization.



HIGH RESISTANCE AND SMALL SIZE

The OP-S8 chair has a completely electrical system for the movement of the footrest and legrest, which is synchronized.

The legrest is concave to hold the patients calves for greater comfort and relaxation.

The range of the backrest movement goes from 80° to the Trendelenburg position.



OPTOMIC®



producing quality

The OPTOMIC OP-S8 chair is designed and manufactured using the best materials and most modern technologies. OPTOMIC's philosophy is based on the idea that both patient and doctor must be able to find a comfortable and safe position, to avoid unnecessary fatigue during an examination or surgical treatment.

CHAIR OP-S8

OP-S8

UP / DOWN MOVEMENT	Electrical
BACKREST AND LEGREST MOVEMENT	Electrical
PROGRAMMABLE MEMORIES	6 (programmable)
MOVEMENT CONTROL	2 keypads on the backrest + Foot-operated joystick
UPHOLSTERY	PUR melting, no seams, latex free, washable
FIRE PROOF CLASS	Fire proof N°UL94HB
FOOTREST	Metallic with automatic fold-down
ARMRESTS	Folding armrests
HEADREST	Variable height and inclination. Removable
	Adjustable headrest optional
CHAIR ROTATION	300° with brake on either side
TRENDELEBURG POSITION	-10° inclination
BACKREST ANGLE	From 90° to 190°
SEATING HEIGHT	Max. 900 mm - Min. 520 mm
MAINS CONNECTION	110/240v - 50-60Hz
DIMENSIONS (BETWEEN BASE TRIMS)	600mm (w) x 730mm (h)
MAXIMUM LOAD	200Kg

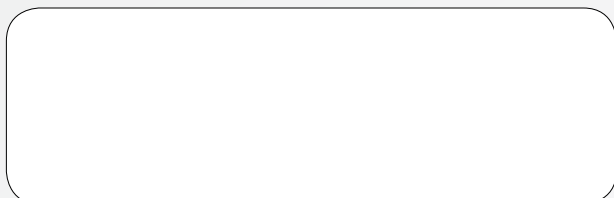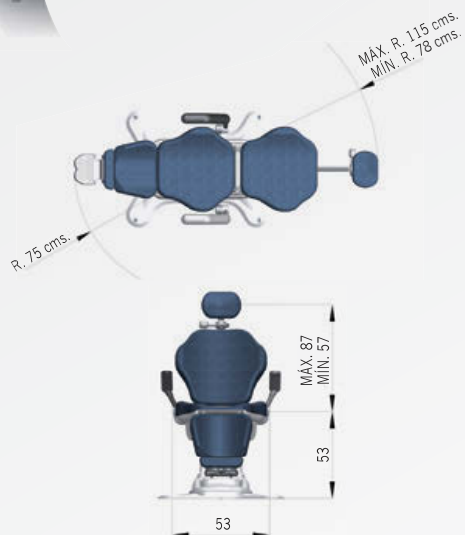

producing quality



Specifications are subject to change without notice

